



Mini-implants/ Mini-implants




B-DENTAL
ORTHODONTIC SUPPLIES

INFINITAS INTRO KIT

DB10-0063 Sterile implants
DB10-0064 Non-Sterile implants

The Infinitas Intro Kit contains:

- Infinitas sterilisation tray
- Standard handle
- Standard inserts: Screwdriver, Mucotome, Bone punch
- 5 mini implants (includes 1 of each non sterile implant)

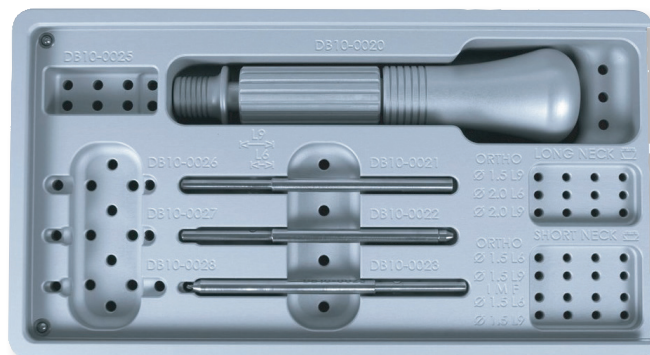


MINI IMPLANT

Blue - Sterile

DB10-0010 1 pc sterile implant

- Diameter - 1.5mm
- Body - 9mm
- Neck - Short
- Typically used for Buccal sites

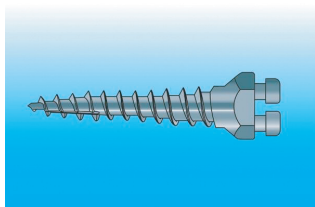


MINI IMPLANT

Green - Sterile

DB10-0016 1 pc sterile implant

- Diameter - 2mm
- Body - 9mm
- Neck - Long
- Typically used for Edentulous areas

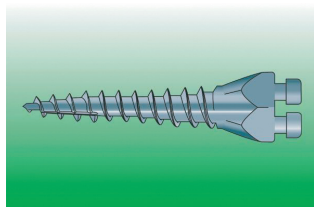


MINI IMPLANT

Magenta - Sterile

DB10-0009 1 pc sterile implant

- Diameter - 1.5mm
- Body Length - 6mm
- Neck - Short
- Used for Incisor sites

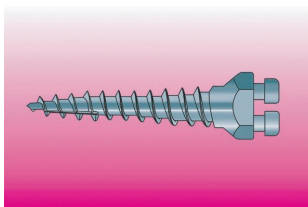


MINI IMPLANT

Purple - Sterile

DB10-0012 1 pc sterile implant

- Diameter - 1.5mm
- Body - 9mm
- Neck - Long
- Typically used for Alveolar Sites

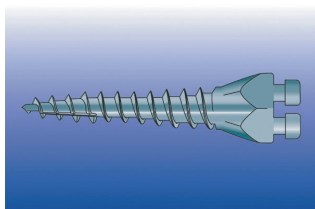


MINI IMPLANT

Yellow - Sterile

DB10-0015 1 pc sterile implant

- Diameter - 2mm
- Body - 6mm
- Neck - Long
- Typically used for Alveolar Sites

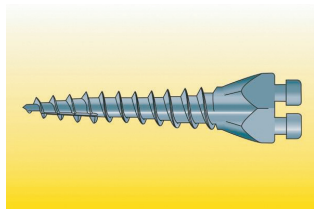


Voor meer informatie gelieve contact op te nemen met B-Dental.

Pour plus d'informations, veuillez contacter B-Dental.

Mini-implants/Mini-implants

INFINITAS MINI IMPLANT

Selecting the correct Infinitas Mini Implant. The Infinitas System can be used to provide bone anchorage in virtually any maxillary/mandibular interproximal site and in the hard palate, although the most frequently utilised insertion site is buccal and mesial to the maxillary first molar. The insertion site is determined by individual anchorage demands and anatomical factors (e.g. bone quantity, root and sinus positions). The exact insertion site is selected by studying each patient’s dental model and radiographs (OPT and/or periapical views) together. General guidelines for usage are:

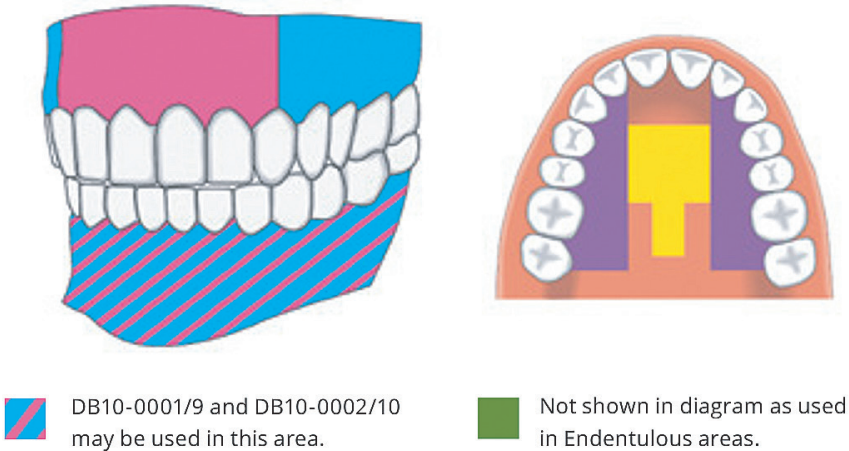

The insertion sites shown are intended for guidance only.

Colour coded	Diameter (mm)	Body Length (mm)	Neck	Non-sterile	Sterile	Typical Insertion sites
	1.5	6	Short	DB10-0001	DB10-0009	Anterior maxilla - perpendicular insertions Mandible - perpendicular insertions
	1.5	9	Short	DB10-0002	DB10-0010	Anterior maxilla - oblique insertions, Posterior maxillary buccal alveolus, Mandible - oblique insertions
	1.5	9	Long	DB10-0004	DB10-0012	Maxilla - palatal alveolus
	2.0	6	Long	DB10-0007	DB10-0015	Mid-palate
	2.0	9	Long	DB10-0008	DB10-0016	Edentulous areas Temporary restorative abutment

